

specifications

The medical gas outlet shall be an Amico Alert-1 series.

Outlets shall be manufactured with a 6 - 1/2" [165 mm] length type "K" 1/2" [12.7 mm] outside diameter (3/8" nominal) size copper inlet pipe stub, which is silver brazed to the outlet body. The body shall be of 1 - 5/16" [33 mm] diameter, one piece brass construction. For positive pressure gas services, the outlet shall be equipped with a primary and secondary check valve and the secondary check valve shall be rated at maximum of 200 psi [1,379 kPa] in the event the primary check valve is removed for maintenance. Outlet bodies shall be gas specific by indexing each gas service to a gas specific dual pin indexing arrangement on the respective identification module.

A large color coded front plate shall be used for ease of gas identification and aesthetic appeal.

A one piece chromed fascia plate shall frame the outlet. With the back rough-in mounted the outlet shall adjust from 1/2" [13 mm] to 1" [25 mm] variation in wall thickness.

The outlets shall be of modular design and include a gas specific 1.6 mm [16 ga.] Steel mounting plate designed to permit on-site ganging of multiple outlets, in any order, on 5" [127 mm] spacing .

The latch-valve assembly shall be D.I.S.S. (Diameter Index Safety System) and shall only accept corresponding gas specific adapters.

All M.R.I. outlets are manufactured from nonferrous materials.

All outlets shall be cleaned and degreased for medical gas service, factory assembled and tested.

The outlets are UL Listed and CSA approved.

Amico products comply with NFPA-99.




features

- ▶ Flow rate exceeds the requirements of NFPA and CSA
- ▶ Accepts only DISS gas specific adapters
- ▶ 100 % hydrostatic testing of each outlet and cleaned for medical gas service
- ▶ Indexed to eliminate interchangeability of gas services
- ▶ Up to 1" [25mm] wall thickness adjustment
- ▶ 360° swivel inlet pipe for easy installation
- ▶ Universal rough-in assembly to accept Quick Disconnects or DISS front adapters (Can be interchanged at any time)
- ▶ Outlets are also available for M.R.I. use

project

wall outlet DISS

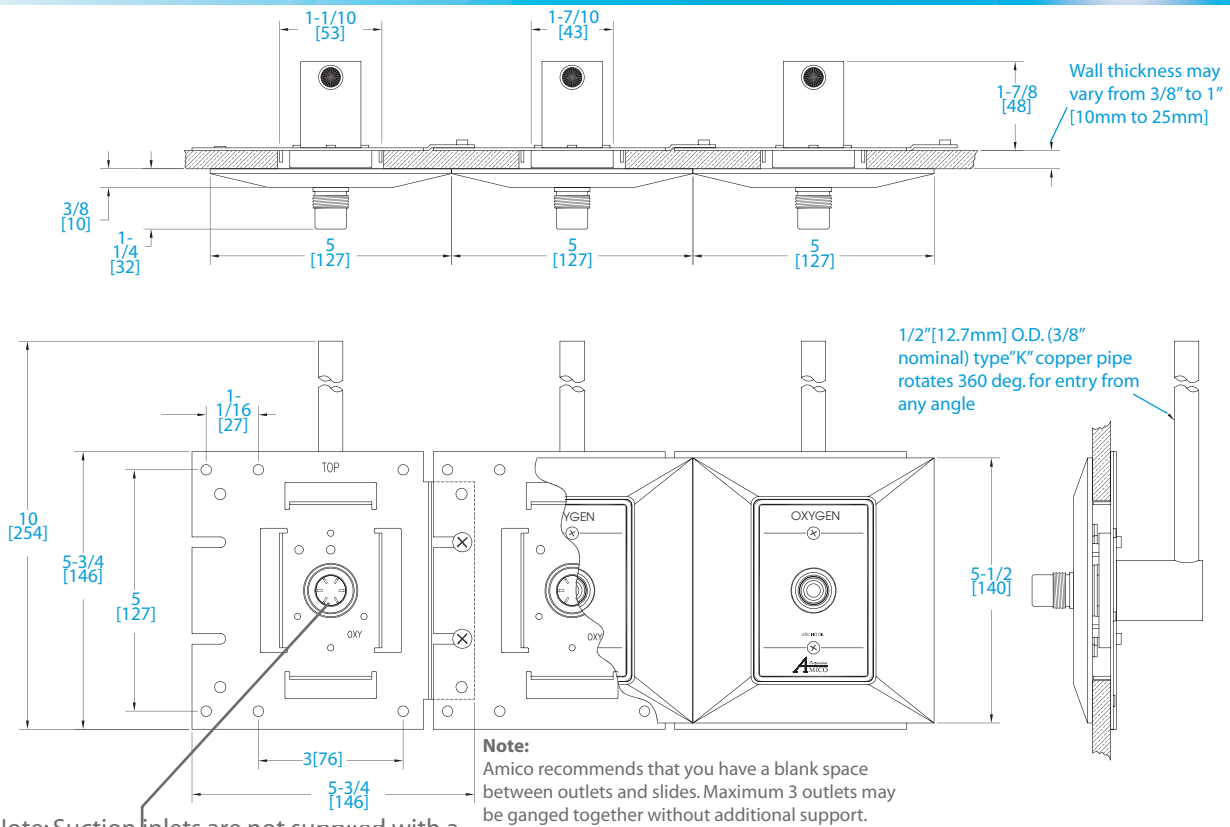

1Note: Suction inlets are not supplied with a secondary check valve. A pressure plug is provided for testing purposes, rated at 150 psi [1,034 kPa]

model numbers

The XXX defines the Gas:

Oxygen	=	OXY
Medical Air	=	AIR
MedVac1	=	VAC
Nitrous Oxide	=	N2O
Nitrogen	=	NIT
Carbon Dioxide	=	CO2
Waste Anesthetic Gas Disposal1 (NFPA)	=	WAG
Anesthetic Gas Scavenging System1 (CSA)	=	AGS

1Please see note above

O-DISWAL-L-XXX

The L defines the Language:

- U = English (NFPA) Oxy, Air, Vac and WAG only
- E = English (CSA/ISO)
- F = French (CSA/ISO)
- S = Spanish (NFPA)

Example:

English NFPA - Oxygen = O-DISWAL-U-OXY

M.R.I. model numbers

O-DISWALM-L-XXX

Example:

English NFPA(M.R.I.) - DISS, Oxygen = O-DISWALM-U-OXY

represented by: