



FAIRFIELD COMPATIBLE RAIL

SPECIFICATIONS:

Each Fairfield requires a minimum of 2 mounting brackets


- 2.0' (609) W-RAIL-F-024
- 3.0' (914) W-RAIL-F-036
- 4.0' (1219) W-RAIL-F-048
- 5.0' (1524) W-RAIL-F-060
- 6.0' (1829) W-RAIL-F-072
- 8.0' (2438) W-RAIL-F-096
- 10.0' (3048) W-RAIL-F-120

The Fairfield Compatible Rail is a complete wall mounted equipment storage system. Installed at each patient location, it provides ready access to commonly used equipment such as suction regulators, flowmeters, exam lights, shelves and collection bottles.

The Fairfield Compatible Rail accommodates a wide variety of snap-in brackets which interface with various equipment. The FCR is permanently attached to the wall in the vicinity of the Medical Gas Outlets. The smooth unobstructed design of the rail's lower surface facilitates easy cleaning.

Ordered separate, each FCR includes aluminum mounting brackets which are field attached to the building structure by the customer. The FCR is available in a variety of lengths and can be field modified to job-site conditions, if required. The FCR is provided with Antique White color laminate inlay, other colors are available upon request.


Extended Wall Mount 1.5"

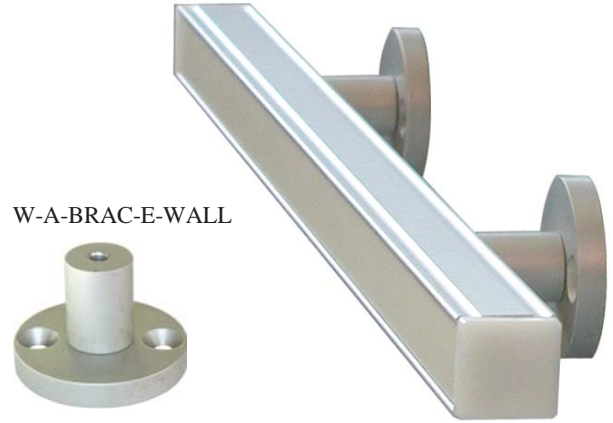

W-A-BRAC-E-WALL



Direct Wall Mount

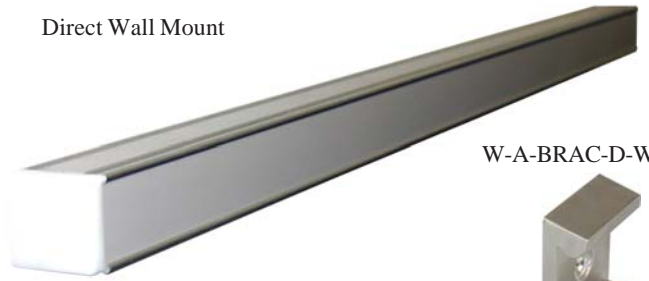

W-A-BRAC-D-WALL



PROJECT:

D I S T R I B U T E D B Y :

A Corporation
AMICO 14 Madison Rd., Fairfield, NJ 07004 USA
85 Fulton Way, Richmond Hill, ON L4B 2N4 Canada

Tel: 1-877-GOAMICO (1-877-462-6426)
Fax: (905) 764-0862 www.amico.com

MAR-03