

SPECIFICATIONS:

The Floor Pedestal shall be an Amico **Alert-1** series.

The floor pedestal shall be made of stainless steel with a #4 satin finish.

It shall be capable of holding medical gas outlets and electrical devices within the structure of the unit.

The electrical devices are provided with electrical boxes and 18" conduit and wire.

There are (2) standard sizes to choose from:

Small Shell - 12" * 10" * 20"



Large Shell - 18" * 10" * 20"

The Floor Pedestal is **UL** Listed to U.S. and Canadian safety standards.



FEATURES:


- Stainless steel shell for aesthetic appearance
- All devices and gas outlets are accessible through the access panel on the back of the pedestal or the removable top
- Custom sizes are available upon request

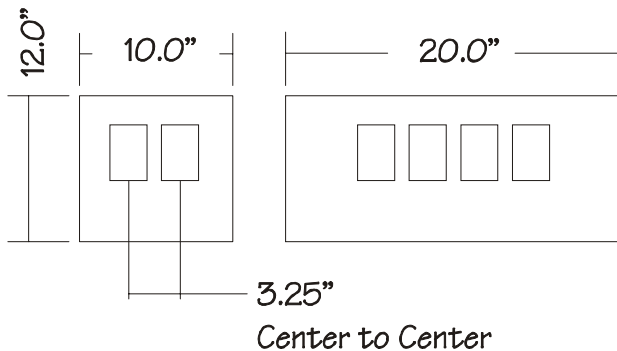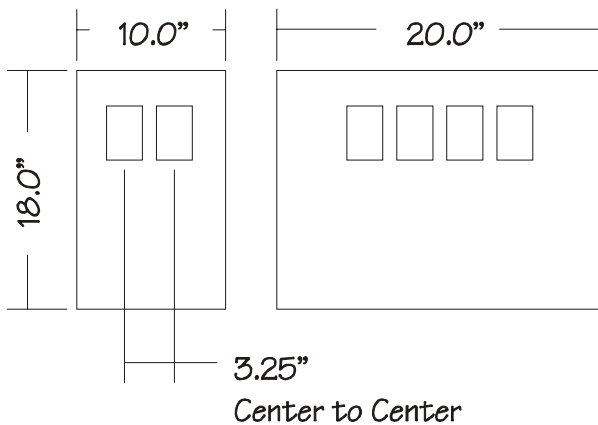


PROJECT:

INSTALLATION INFORMATION

Small Pedestal

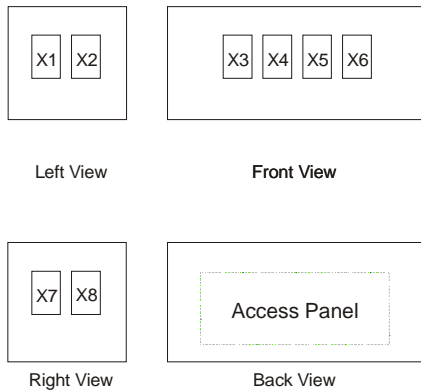

Large Pedestal

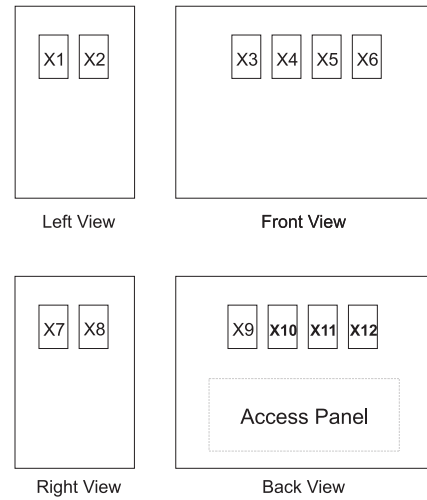

EQUIPMENT LAYOUT

X1	
X2	
X3	
X4	
X5	
X6	
X7	
X8	
X9	
X10	
X11	
X12	

INSTALLATION:

Use (four) 3/8"-1/2" diameter mounting bolts (provided by contractor) to mount the pedestal securely on the finished floor. Remove the top cover and access panel for piping and wiring.

MODEL NUMBER

Small: W-PEDESTAL-12

Large: W-PEDESTAL-18

What style of outlet do you require:
DISS, Ohmeda, Chemetron or Puritan
Bennett

FLOOR PATTERN DIMENSIONS:

