



Installation Instructions

noble series
neonatal intensive care unit (NICU)




On-site Connections

Pipe connections between base and top, to be made using fittings and pipes provided by Amico.

All conduits from the base to be carried through the chase to the top of the unit.

Monitor to be connected directly to the monitor back box, through the chase.



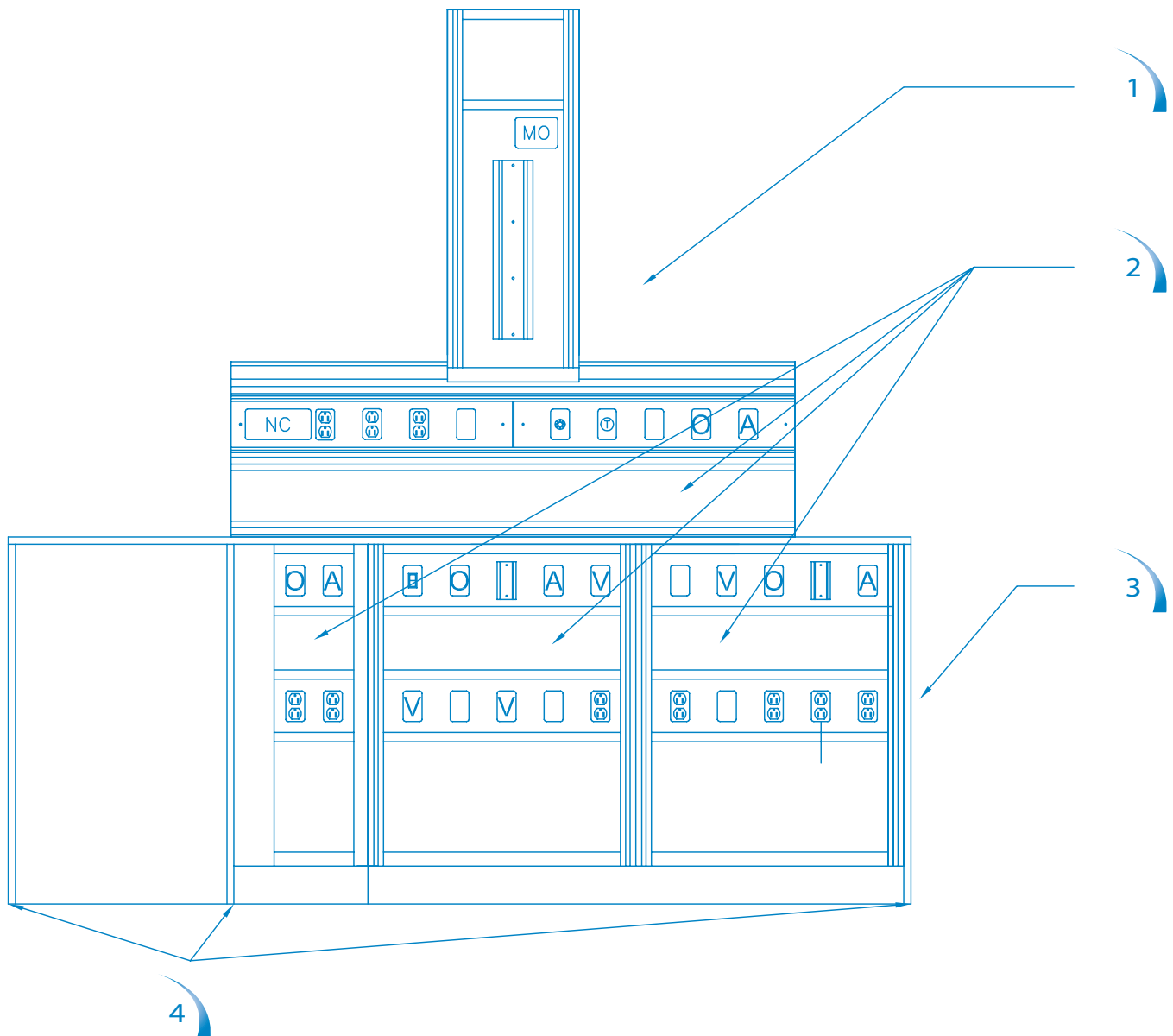

1 Attach top unit to the base using holes (4 holes) and hardware provided.

2 Remove Formica panels for access to mounting holes.

3 Base assembly is pre-piped and pre-wired.

4 Remove white plastic feet for permanent installation only.



85 Fulton Way, Richmond Hill
Ontario, L4B 2N4, Canada

14 Madison Road, Fairfield
NJ 07004, USA

Toll Free Tel: 1.877.462.6426
Toll Free Fax: 1.866.440.4986
Tel: 905.764.0800 | Fax: 905.764.0862
www.amico.com

